

Mini Masaj Aleti

Ürün Tanımı

Bu cihaz; geleneksel Çin tıp prensiplerine göre, uluslararası gelişmiş mikrobilgisayar elektronik teknolojisi kullanılarak yüksek teknoloji ile üretilmiştir. Geleneksel Çin tıbbi masaj tekniklerini simüle eden bir mikrobilgisayar ile çalışır.

Masaj yapar, yaşam enerjisi ve kanı dengeler, yorgunluğu giderir, vücudu güçlendirir


Ürün Özellikleri

- Çeşitli masaj modları arasında seçim yapılmasına imkan sunar.
- Omuzlar, bel, dizler, eller, ayaklar ve bacaklarda kullanıma uygundur.
- Kullanım kolaylığı sağlar.
- Herhangi bir yan etki veya güvenlik tehdidi olmadan tedavi sağlar.

1

Mini Masaj Aleti

Ürün Yapısı

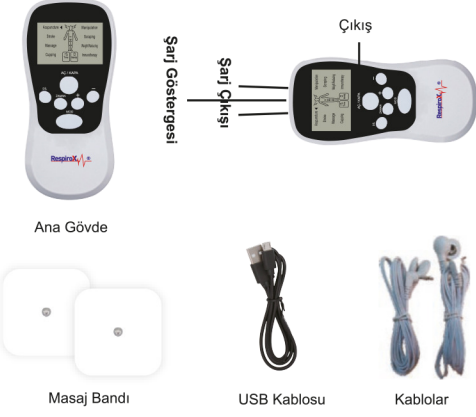

Uyarı

- Sağlığınız için kullanıma başlamadan önce mutlaka doktorunuza danışın.
- Aşağıda belirtilen kişilerde kullanıma uygun değildir:
 - Kötü huylu tümör hastaları
 - Hamile veya regl kadınlar
 - Kalp krizi geçirmiş kişiler
 - Akut hastalığı olan kişiler
 - Doktor kontrolünde olan hastalar (fiziksel olarak rahatsız hissediyorlarsa)

Not: Bu kılavuzda gösterilen aksesuarlar referans amaçlıdır. Ürün ile birlikte sağlanacak aksesuarlar kullanıcı tarafından belirlenir.

2

Mini Masaj Aleti

Kullanım Şekli

1. Kablunun bir ucunu cihazdaki uygun çıkışa diğer ucunu ise masaj bandı düzgün bir şekilde takın.
2. Masaj bandının tedavinin uygulanacağı alana yapıştırın. (yapıştırmadan önce alanın temiz olmasına dikkat edin)
3. Cihazı çalıştırmak için cihazın sol üst tarafında bulunan açma tuşuna basın.
4. Cihazda sekiz farklı masaj modu bulunmaktadır.
5. Cihaz çalışırken masaj levhalarının her iki kutbunu da düzgün bir şekilde yerleştirin. Levhalar tam olarak veya hiç yapışmazsa cihaz otomatik olarak kapanacaktır.
6. Bu cihazda 5-15 dakika otomatik kapanma özelliği bulunmaktadır. LED ekranda bulunan "time" alanında zaman zaman kapanmaya kalan süre görüntülenir. "0-time" görüntülendiğinde cihaz otomatik olarak çalışmayı durdurur.

Hata Analizi ve Sorun Giderme

Hata	Sebebi	Çözüm
Otomatik Kapanma	Bant tam yapışmamış veya kablo tam takılı değil	Masaj yapılan bölgeni temizleyip bantı tekrar takın, kablounun doğru şekilde takılı olduğunu kontrol edin.
Fiziksel terapi etkisi zayıf	Çalışma hatası ve voltajda dalgalanma	Cihazı kapatıp tekrar çalıştırın
Vücudun bazı bölgelerinde rahatsızlık hissi	Çok uzun süreli kullanım	Tek seferde 30 dakikadan fazla kullanmayın
	Aynı yerde çok uzun süre kullanım	Kullanım bölgesini değiştirin

3

Mini Masaj Aleti

Notlar

- Kullanım öncesinde masaj bantlarını yerleştirin ardından cihazı çalıştırın.
- Kullanım sırasında bantların yerini değiştirmeyin.
- Kullanım sırasında rahatsızlık hissederseniz cihazı kapatın.
- Şebeke gerilimi dışında bir güç kaynağı kullanmayın.
- Cihazı tek başınıza onarmaya ve modifiye etmeye çalışmayın.
- Banyo ve mutfak benzeri ıslak ve/veya sıcak yerlerde kullanmayın.

Depolama

- Cihazda kir veya toz varsa nötr deterjanla batırılmış yumuşak havlu ile temizleyin. Havalanarak kendi kendine kurumasını sağlayın.
- Cihaz uzun süre kullanılmayacaksa pilini çıkararak muhafaza edin.
- Cihazı sıcak vb. sıcak yüzeylerde muhafaza etmeyin.
- Cihazı nemli ortamlarda muhafaza etmeyin.



UYARI

Cihazı temizlemek için sülfürik asit vb. aşındırıcı sıvılar kullanmayın.

4

Mini Masaj Aleti

Satış Sonrası Hizmetler ve Üretici Sorumlulukları

1. Ürün garanti süresi satın alma tarihinden itibaren bir yıldır.
2. Garanti süresi boyunca oluşacak arızalarda, kullanıcı kusuru olmamak şartıyla, fatura ve garanti belgesi ile ücretsiz bakım ve onarım talep edilebilir. Aşağıdaki durumlarda, garanti süresi devam eden ürünler dahil olmak üzere, malzeme ve servis ücreti talep edilebilir:
 - Fatura ve garanti belgesi yoksa
 - Yanlış kullanımdan kaynaklı hatalarda
 - Kullanıcı müdahalesi olan durumlarda
 - Mücbir sebepler ile oluşan hasarlarda

1. Ürün ile ilgili tüm testler gerçekleştirilmiş ve kullanımının güvenli ve uygun olduğu kanıtlanmıştır.
2. Cihazı doğru bir şekilde kullanabilmek için bu kılavuzu dikkatli bir şekilde okuyun ve saklayın.
3. Firma ürün ve özelliklerini önceden bildirim yapmadan değiştirme ve geliştirme hakkını saklı tutar.

Adaptör ile Kullanım



Pil ile Kullanım



5

3 adet AAA pil



Shenzhen AXD Electronic Co., Ltd.

Building A2, Fuyong Zhihui Garden, Baoan District, Shenzhen City, Guangdong Province

İTHALATÇI: Ece Tıbbi Cihazlar ve Medikal San. Tic. A.Ş.

Maltepe Mah. Hastane Yolu Sok. No: 6 Zeytinburnu, İstanbul, Türkiye

www.respirox.com.tr / www.ecemedikal.net

Tel: 444 78 96

Yayın Tarihi: 09.12.2022

6



Çok Fonksiyonlu Çift Girişli Masaj Aleti

RTC-01

- Masaj yapar
- Enerjisi dengeler
- Zinde hissettirir
- Yorgunluğu giderir

Kullanım Kılavuzu



Kurulum ve Kullanım Videosu için For installation and Usage Video Taratınız. Scan it.



Mini massager

Product description

According to the traditional Chinese medicine principle, using the international advanced microcomputer electronic technology, the research and development center developed a kind of multifunctional microcomputer electronic massage product with high-tech content. It is controlled by micro-computer under the low-cycle simulation of traditional Chinese medicine massage techniques, through the human body Meridian system conditioning the body, play through the channels and collaterals to harmonize Qi and blood, Yin into coordination, health care role.

Dredge channels and collaterals, harmonize Qi and blood, relieve fatigue and strengthen body

Product features

- * Various massage modes, such as stroke, intelligence, massage, cupping, massage, scrapping, slimming, immunomodulation, etc.
- * Suitable for shoulders, waist, knees, hands and feet. Legs and everything.
- * Easy to use, easy to operate, self-use, suitable for giving away.
- * Purely natural therapy, with no side effects or safety concerns.

1

Mini massager

Product Construction



Fair warning

1. For your health and that of your family, please consult your physician before using.
2. Malignant Tumor Patient
3. Pregnant or menstruating.
4. Heart attack victim
5. Patients with acute disease.
6. People who are under the care of a physician, especially if they feel physically abnormal

2

Mini massager

Method of use

1. Please insert the lead plug into the output Jack of the massager and connect the other end to the massage patch correctly.
2. Please paste the patch on the part that needs massaging (please clean the part that needs massaging first).
3. Long press the on button (upper left side of the product) to open the massager.
4. The machine is equipped with eight massage modes, such as hammer-stroke, intelligence, massage, cupping, massage, scrapping, slimming, immunomodulation, etc.
5. When the product is working, please stick the two poles of the massage placard firmly, in case it is detected that the load is too light or there is no load and the machine will shut down automatically.
6. This machine is equipped with 5-15 minutes automatic shutdown function. The LED screen "time" randomly displays the remaining working time. When the LED screen displays "0-time", the machine automatically stops working.

Failure Analysis and treatment

An unusual situation	The reason	Handle it
Automatic shutdown	The patch is falling off or the connecting wire is not in good contact	Please clean the massage area and reattach the patch to the check patch connector
The effect of physical therapy is weak	Operation Error or voltage instability	Turn off the power and reboot
Some parts of the body feel uncomfortable	If it has been used too long Whether massages the same place for a long time	Use for more than 30 minutes at a time Please replace the other part

3

Mini massager

Points to note

- * Before using, please paste the massage patch on the skin, and then turn on the power.
- * Do not place the patch on or transfer it to another person during use.
- * Stop using the product if you feel unwell during use. * Do not use power supply other than rated voltage
- * Do not disassemble, repair or modify this machine except for our designated repair technician. * Do not use in wet or hot places such as bathrooms and kitchens.

Storage method

- * If there is dirt or dust, please wipe with a soft cloth dipped in neutral detergent and let it dry naturally with ventilation.
- * If the machine has been out of use for too long, please take out the battery and save it separately.
- * Do not store in a hot place such as a stove or kitchen when not in use.
- * Do not store in a damp place when not in use.



Warning

Please do not use the water, sulfuric acid and other corrosive liquid wipe.

4

Mini massager

After-sales service and producer responsibility

1. The quality guarantee period of this product is one year, counting from the date of purchase.
2. During the warranty period, if not damaged by manpower, the shop can enjoy free repair and service with invoice and warranty card. In the following cases, though within the warranty period, material cost or service fee will be charged according to the circumstances:
 - failure to produce purchase invoice and warranty card
 - damage due to misuse or error
 - failure due to self-repair or modification
 - damage due to force majeure

Charging charge



Dry-cell models



Use 3 # 7 dry cells

5



Shenzhen AXD Electronic Co., Ltd.

Building A2, Fuyong Zhihui Garden, Baoan District, Shenzhen City, Guangdong Province

IMPORTER: Ece Tıbbi Cihazlar ve Medikal San. Tic. A.Ş.

Maltepe Mah. Hastane Yolu Sok. No: 6 Zeytinburnu, İstanbul, Türkiye

www.respirox.com.tr / www.ecemedikal.net

Tel: 444 78 96

Date of Issue: 09.12.2022

6

Mini massager

Product description

According to the traditional Chinese medicine principle, using the international advanced microcomputer electronic technology, the research and development center developed a kind of multifunctional microcomputer electronic massage product with high-tech content. It is controlled by micro-computer under the low-cycle simulation of traditional Chinese medicine massage techniques, through the human body Meridian system conditioning the body, play through the channels and collaterals to harmonize Qi and blood, Yin into coordination, health care role.

Dredge channels and collaterals, harmonize Qi and blood, relieve fatigue and strengthen body

Product features

- * Various massage modes, such as stroke, intelligence, massage, cupping, massage, scrapping, slimming, immunomodulation, etc.
- * Suitable for shoulders, waist, knees, hands and feet. Legs and everything.
- * Easy to use, easy to operate, self-use, suitable for giving away.
- * Purely natural therapy, with no side effects or safety concerns.

1

Mini massager

Product Construction



Fair warning

1. For your health and that of your family, please consult your physician before using.
2. Malignant Tumor Patient
3. Pregnant or menstruating.
4. Heart attack victim
5. Patients with acute disease.
6. People who are under the care of a physician, especially if they feel physically abnormal

2

Mini massager

Method of use

1. Please insert the lead plug into the output Jack of the massager and connect the other end to the massage patch correctly.
2. Please paste the patch on the part that needs massaging (please clean the part that needs massaging first).
3. Long press the on button (upper left side of the product) to open the massager.
4. The machine is equipped with eight massage modes, such as hammer-stroke, intelligence, massage, cupping, massage, scrapping, slimming, immunomodulation, etc.
5. When the product is working, please stick the two poles of the massage placard firmly, in case it is detected that the load is too light or there is no load and the machine will shut down automatically.
6. This machine is equipped with 5-15 minutes automatic shutdown function. The LED screen "time" randomly displays the remaining working time. When the LED screen displays "0-time", the machine automatically stops working.

Failure Analysis and treatment

An unusual situation	The reason	Handle it
Automatic shutdown	The patch is falling off or the connecting wire is not in good contact	Please clean the massage area and reattach the patch to the check patch connector
The effect of physical therapy is weak	Operation Error or voltage instability	Turn off the power and reboot
Some parts of the body feel uncomfortable	If it has been used too long Whether massages the same place for a long time	Use for more than 30 minutes at a time Please replace the other part

3

Mini massager

Points to note

- * Before using, please paste the massage patch on the skin, and then turn on the power.
- * Do not place the patch on or transfer it to another person during use.
- * Stop using the product if you feel unwell during use. * Do not use power supply other than rated voltage
- * Do not disassemble, repair or modify this machine except for our designated repair technician. * Do not use in wet or hot places such as bathrooms and kitchens.

Storage method

- * If there is dirt or dust, please wipe with a soft cloth dipped in neutral detergent and let it dry naturally with ventilation.
- * If the machine has been out of use for too long, please take out the battery and save it separately.
- * Do not store in a hot place such as a stove or kitchen when not in use.
- * Do not store in a damp place when not in use.



Warning

Please do not use the water, sulfuric acid and other corrosive liquid wipe.

4

Mini massager

After-sales service and producer responsibility

1. The quality guarantee period of this product is one year, counting from the date of purchase.
2. During the warranty period, if not damaged by manpower, the shop can enjoy free repair and service with invoice and warranty card. In the following cases, though within the warranty period, material cost or service fee will be charged according to the circumstances:
 - failure to produce purchase invoice and warranty card
 - damage due to misuse or error
 - failure due to self-repair or modification
 - damage due to force majeure

Charging charge



Dry-cell models



Use 3 # 7 dry cells

5



Shenzhen AXD Electronic Co., Ltd.

Building A2, Fuyong Zhihui Garden, Baoan District, Shenzhen City, Guangdong Province

IMPORTER: Ece Tıbbi Cihazlar ve Medikal San. Tic. A.Ş.

Maltepe Mah. Hastane Yolu Sok. No: 6 Zeytinburnu, İstanbul, Türkiye

www.respirox.com.tr / www.ecemedikal.net

Tel: 444 78 96

Date of Issue: 09.12.2022

6

Mini massager

Product description

According to the traditional Chinese medicine principle, using the international advanced microcomputer electronic technology, the research and development center developed a kind of multifunctional microcomputer electronic massage product with high-tech content. It is controlled by micro-computer under the low-cycle simulation of traditional Chinese medicine massage techniques, through the human body Meridian system conditioning the body, play through the channels and collaterals to harmonize Qi and blood, Yin into coordination, health care role.

Dredge channels and collaterals, harmonize Qi and blood, relieve fatigue and strengthen body

Product features

- * Various massage modes, such as stroke, intelligence, massage, cupping, massage, scrapping, slimming, immunomodulation, etc.
- * Suitable for shoulders, waist, knees, hands and feet. Legs and everything.
- * Easy to use, easy to operate, self-use, suitable for giving away.
- * Purely natural therapy, with no side effects or safety concerns.

1

Mini massager

Product Construction



Fair warning

1. For your health and that of your family, please consult your physician before using.
2. Malignant Tumor Patient
3. Pregnant or menstruating.
4. Heart attack victim
5. Patients with acute disease.
6. People who are under the care of a physician, especially if they feel physically abnormal

2

Mini massager

Method of use

1. Please insert the lead plug into the output Jack of the massager and connect the other end to the massage patch correctly.
2. Please paste the patch on the part that needs massaging (please clean the part that needs massaging first).
3. Long press the on button (upper left side of the product) to open the massager.
4. The machine is equipped with eight massage modes, such as hammer-stroke, intelligence, massage, cupping, massage, scrapping, slimming, immunomodulation, etc.
5. When the product is working, please stick the two poles of the massage placard firmly, in case it is detected that the load is too light or there is no load and the machine will shut down automatically.
6. This machine is equipped with 5-15 minutes automatic shutdown function. The LED screen "time" randomly displays the remaining working time. When the LED screen displays "0-time", the machine automatically stops working.

Failure Analysis and treatment

An unusual situation	The reason	Handle it
Automatic shutdown	The patch is falling off or the connecting wire is not in good contact	Please clean the massage area and reattach the patch to the check patch connector
The effect of physical therapy is weak	Operation Error or voltage instability	Turn off the power and reboot
Some parts of the body feel uncomfortable	If it has been used too long Whether massages the same place for a long time	Use for more than 30 minutes at a time Please replace the other part

3

Mini massager

Points to note

- * Before using, please paste the massage patch on the skin, and then turn on the power.
- * Do not place the patch on or transfer it to another person during use.
- * Stop using the product if you feel unwell during use. * Do not use power supply other than rated voltage
- * Do not disassemble, repair or modify this machine except for our designated repair technician. * Do not use in wet or hot places such as bathrooms and kitchens.

Storage method

- * If there is dirt or dust, please wipe with a soft cloth dipped in neutral detergent and let it dry naturally with ventilation.
- * If the machine has been out of use for too long, please take out the battery and save it separately.
- * Do not store in a hot place such as a stove or kitchen when not in use.
- * Do not store in a damp place when not in use.



Warning

Please do not use the water, sulfuric acid and other corrosive liquid wipe.

4

Mini massager

After-sales service and producer responsibility

1. The quality guarantee period of this product is one year, counting from the date of purchase.
2. During the warranty period, if not damaged by manpower, the shop can enjoy free repair and service with invoice and warranty card. In the following cases, though within the warranty period, material cost or service fee will be charged according to the circumstances:
 - failure to produce purchase invoice and warranty card
 - damage due to misuse or error
 - failure due to self-repair or modification
 - damage due to force majeure

Charging charge



Dry-cell models



Use 3 # 7 dry cells

5



Shenzhen AXD Electronic Co., Ltd.

Building A2, Fuyong Zhihui Garden, Baoan District, Shenzhen City, Guangdong Province

IMPORTER: Ece Tıbbi Cihazlar ve Medikal San. Tic. A.Ş.

Maltepe Mah. Hastane Yolu Sok. No: 6 Zeytinburnu, İstanbul, Türkiye

www.respirox.com.tr / www.ecemedikal.net

Tel: 444 78 96

Date of Issue: 09.12.2022